



THE ORTHODONTIST'S LIGHT

AO
AMERICAN
ORTHODONTICS



BLUE RAY 3 THE ORIGINAL ORTHODONTIST'S CURING LIGHT

The Blue Ray 3 curing light is powerful, customizable, and made specifically for the orthodontist, delivering efficient and dependable curing every time.

POWERFUL CURING

Blue Ray 3 delivers up to 3000 mW/cm² of curing power and is one of the most powerful lights in the industry.


ORTHO FRIENDLY TIMING MODES

Customizable timing options allow one-button curing of a few brackets or an entire arch. Choose 3+3 second mode for high speed curing or 5+5 second mode for extra bonding security.

3-YEAR WARRANTY

Cure with confidence knowing that your Blue Ray 3 is covered by a 3-year warranty on the hand piece, battery, and charging base.

- Available in brushed stainless steel or polished black finish
- Four customizable alert cues, including visual microflash
- Ergonomic, lightweight design for ease of use
- One charge delivers 45 minutes of continuous light, enough for 12-15 patients
- LED light won't lose power over time and emits no UV wavelengths
- Hygienic sealed case
- Power verification window on charging base ensures proper curing power
- Kit includes both molar-shaped and oval eye shields

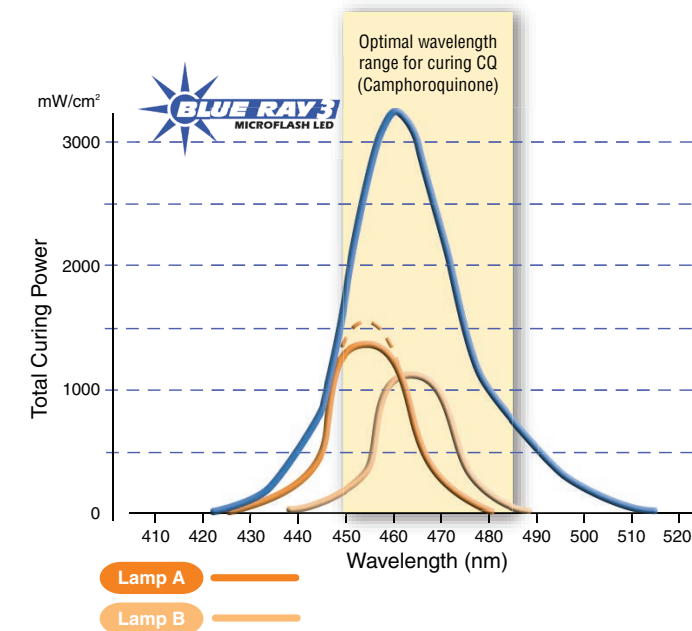


KIT

Includes: hand piece, battery, molar-shaped eye shield, oval eye shield, charging base, power cord, and 5.5 mm light guide


Black Stainless Steel kit	851-945
Silver Stainless Steel kit	851-940

REORDER COMPONENTS

(1) Hand piece - Black	851-946
Hand piece - Silver	851-941
(2) Battery - Black	851-947
Battery - Silver	851-932
(3) Charging base/cradle - Black	851-948
Charging base/cradle - Silver	851-935
(4) Molar-shaped eye shield	851-933
(5) Oval eye shield	851-934
Power Cord	851-936
(6) 5.5 mm light guide - black	851-913B
5.5 mm light guide - frosted (optional)	851-913F
7.5 mm light guide - frosted (optional)	851-915F

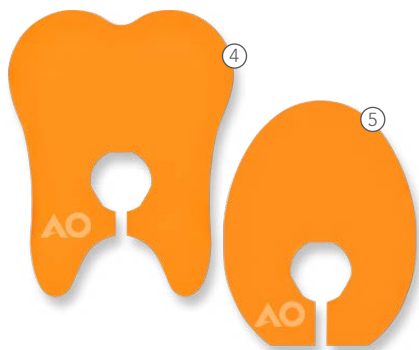

Close-up
of Hand Piece

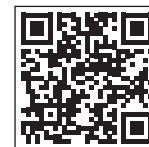

Scan for more
information.



3524 Washington Avenue
Sheboygan, Wisconsin, 53081 USA

Phone: +1 920 457 5051
Email: info@americanortho.com
Web: americanortho.com

©2015 American Orthodontics Corporation
All Rights Reserved

FOLLOW US ON



3 SECOND MODE

CERAMIC BRACKETS	3 seconds through face of bracket
METAL BRACKETS	3 seconds mesial and 3 seconds distal
DIRECT BOND TUBES	3 seconds mesial and 6 seconds occlusal
MOLAR BANDS	3 seconds per cusp = 12 seconds total

5 SECOND MODE

CERAMIC BRACKETS	5 seconds through face of bracket
METAL BRACKETS	5 seconds mesial and 5 seconds distal
DIRECT BOND TUBES	5 seconds mesial and 10 seconds occlusal
MOLAR BANDS	5 seconds per cusp = 20 seconds total